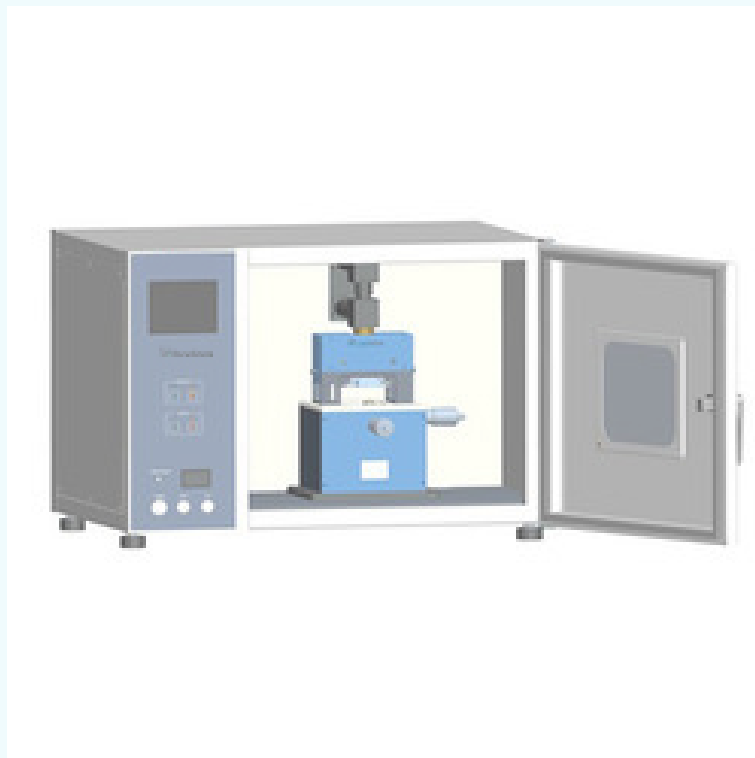




Atomic force microscope LD-LAFM-A11

Labodam Equipment Ltd.
info@labodam.com | www.labodam.com

Atomic force microscope LD-LAFM-A11

Overview

Atomic force microscope LAFM-A11 comes with combined design of scan head and sample stage, to give strong anti-vibration performance. Adopted with servomotor, drives the sample approaching tip manually or automatically, to evaluate precise scanning area position. Equipped with precision laser detection and probe alignment device, make easy adjustment of laser beam.

Features:

- Large sample transfer range
- Optical observation system for checking tip & sample's position
- Modular electronic system for easy maintenance
- CCD observing system
- Equipped with servomotor to achieve CCD auto focusing
- Provides highly accurate results

Applications:

It finds best solution for applications in biochemistry for tissues, cells, cellular components imaging, nanotechnology for imaging of polymers, nanomaterials, and in chemistry physics for imaging of surface metals elements.

Specifications:

Scan rate	0.6 Hz to 4.34 Hz
Net weight	50 kg
Resolution	X/Y axis: 0.2 nm, Z axis: 0.05 nm
Scan angle	0 to 360°
Sample size	Ø? 90 mm, H? 20 mm
Gross Weight	87.4 kg
Feedback type	DSP digital feedback
Operation modes	Contact mode, friction mode, extended modes of tapping phase, MFM, EFM, scan angle 0 to 360°
PC connections:	USB 2.0
Sample movement	0 to 20 mm
Scanning control	XY: 18 bit D/A & double 16 bit A/D multiple channel simultaneously
Windows software	Compatible with windows 98/2000/XP/7/8

Maximum scan range	X/Y axis: 50 μm , Z axis: 5 μm
Instrument Diemnsion	700 $\times$ 500 $\times$ 460 mm
Feedback sampling rate	64 KHz
Types of sampling pixel	256 $\times$ 256, 512 $\times$ 512
Pulse width of approaching motor	10 $\pm$ 2 ms
Optical system/ Magnification of CCD	Magnification: 10x, Resolution: 1 μm